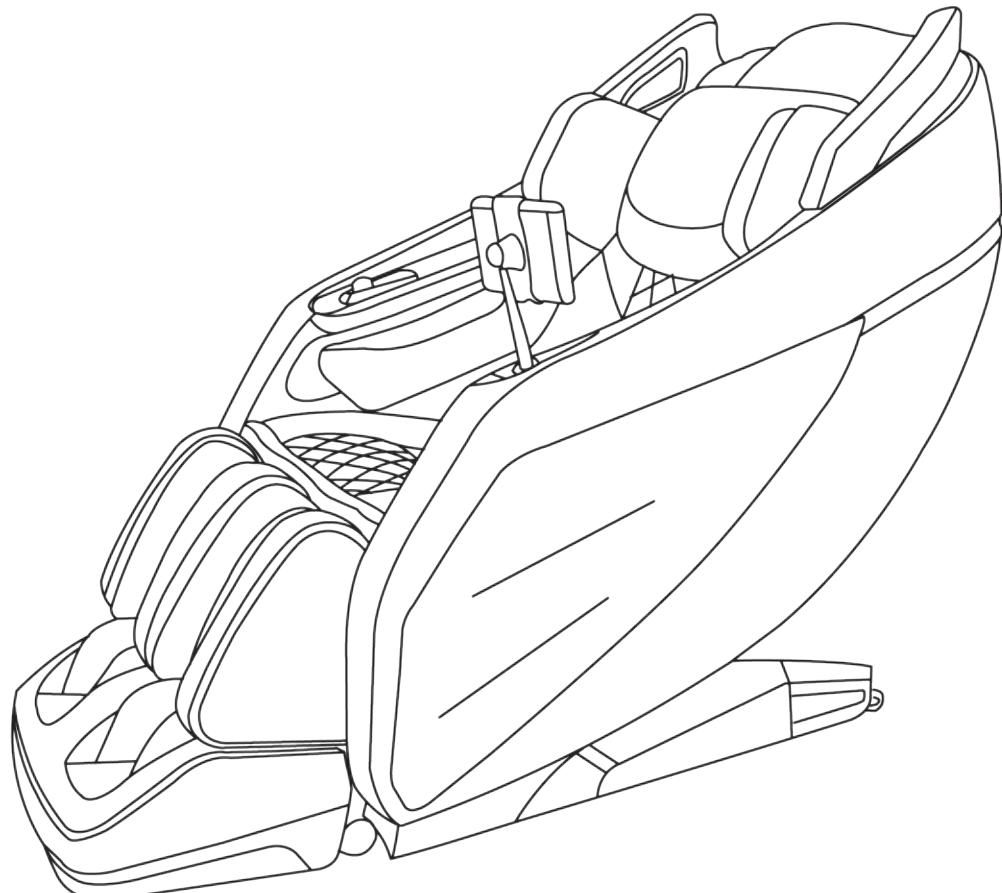


Osaki

OS-3D Champ 2 INSTRUCTION MANUAL

Enjoy the space capsule Massage Chair



PLEASE BE SURE TO READ THIS MANUAL CAREFULLY TO LEARN THE OPERATIONS OF THIS PRODUCT. DO NOT DISPOSE OF THE MANUAL AND KEEP FOR FUTURE REFERENCE.

WARRANTY SERVICE

Email: Service@OsakiTitan.com

Phone: 1-888-848-2630 ext.3

SCAN TO SUBMIT TICKET:



You can also submit ticket directly using below link:

<https://titanchair.freshdesk.com/support/home>

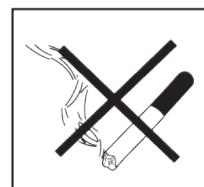
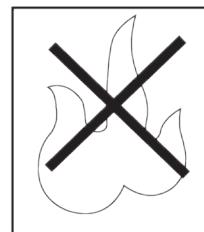
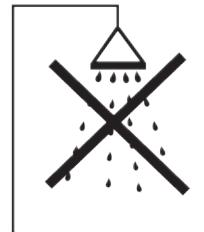
CONTENTS

IMPORTANT SAFETY INSTRUCTIONS	1-4
FUNCTION OF EACH PART.....	5
OPERATION.....	6-20
SAFETY INSTRUCTIONS.....	21
DAILY CARE AND MAINTENANCE	21
TROUBLE SHOOTING	22
SPECIFICATION.....	22

IMPORTANT SAFETY INSTRUCTIONS

OPERATING SURROUNDINGS

- Do not use or store the massage chair in a damp environment : i.e. bathroom, sauna, spa to avoid any mold, electric shock, fire or mechanical failure.
- Do not use the massage chair outdoors
- Do not use the massage chair under dusty, smeary, or anoxic conditions
- Do not place the massage chair near a heater, stove, or near sunlight radiation to avoid inner components overheating and causing fire
- Make sure to place the massage chair on a flat level ground
- The bracket part should be placed on the plain ground
- To avoid electric shock, please make sure to connect an earth wire and have it grounded.
- Don't use the machine while smoking



FOLLOWING CONDITIONS SHOULD CONSULT DOCTORS BEFORE USE

- Suffered from malignant tumor
- Serious heart disease
- Serious skin disease
- Serious boned hyperplasia
- During Pregnancy
- If there is a bend in the back caused by an accident or disease
- Arthritis Disease
- Implanted pacemaker for medical purpose
- Suffered from bowels disease (i.e. gastritis, entritis, diarrhea, hepatitis)
- Serious blood circulation trouble
- Consciousness trouble
- Frequent high fever
- Please consult your doctor for any reason prior to use for the best medical advise

⚠ WARNING

Children under the age of 18 or people with disabilities (physical and mental) are not allowed to operate the massage chair without a doctors consent to prevent cases of electric shock, injuries, or mechanism malfunctions.

IMPORTANT SAFETY INSTRUCTIONS

For safety reasons and in order to clearly indicate potential dangers should the product be used improperly, instructions have been marked in the following manners and should be strictly noted

🚫 CAUTIONS

Non-professional technicians are not allowed to dismantle or repair the chair to avoid fire, electrical shock or serious harm



Do not spray or directly pour liquids to the switch when cleaning to avoid insulation fault, electrical shock or short-circuiting.



When moving the chair in the process of storing or maintenance, be sure to disconnect the power plug from socket



Do not detach the plug with wet hands

Make sure to shut off and detach the power plug when cleaning or storing product



During cleaning or moving the machine, make sure to detach power cord to avoid fire or electrical shock

Do not place the machine onto an unstable surface to avoid roasting, insulation failure, electrical shock or turnover



When storing the machine, make sure to disconnect power plug to avoid fire caused by electricity leaking

Clean power plug periodically to avoid fire due to dust



Do not use it in any unsafe places or for other purposes



It is the main reason of product damage or accident

Do not store the machine onto an unstable surface



Unexpected earthquake or external force will cause damage or injury

Keep supply cord unobstructed at all time, i.e. not twisted, knotted, over flexed or pulled and do not place anything heavy on it



The appliance has a heated surface. Persons sensitive to heat must be careful when using the appliance.



IMPORTANT SAFETY INSTRUCTIONS

ⓧ CAUTIONS

Do not modify the structure or change components by unauthorized person

Do not dismantle the machine by unauthorized person to avoid electrical shock, short-circuit or fire



Do not store it near a fire or flammable substances

It may cause product manipulation, deviation, or damage.



Should any abnormality occur(i.e. scorching smell), stop using immediately and detach power plug from socket

Operating under an abnormal situation will cause fire or electrical shock, please contact your local seller or authorized service agent for maintenance



Only rated voltage must be used



Stop using immediately if supply cord or power plug be damaged or loosen



Be sure to keep children away during operating to avoid any hazard



Keep away from direct sunlight



If the supply cord is damaged, it must be replaced by the manufacturer or its authorized persons for repair or maintenance to avoid any hazard



IMPORTANT SAFETY INSTRUCTIONS

RF exposure statement

This equipment complies with the FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

FCC Warning

This device complies with Part 18 and Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE 1: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

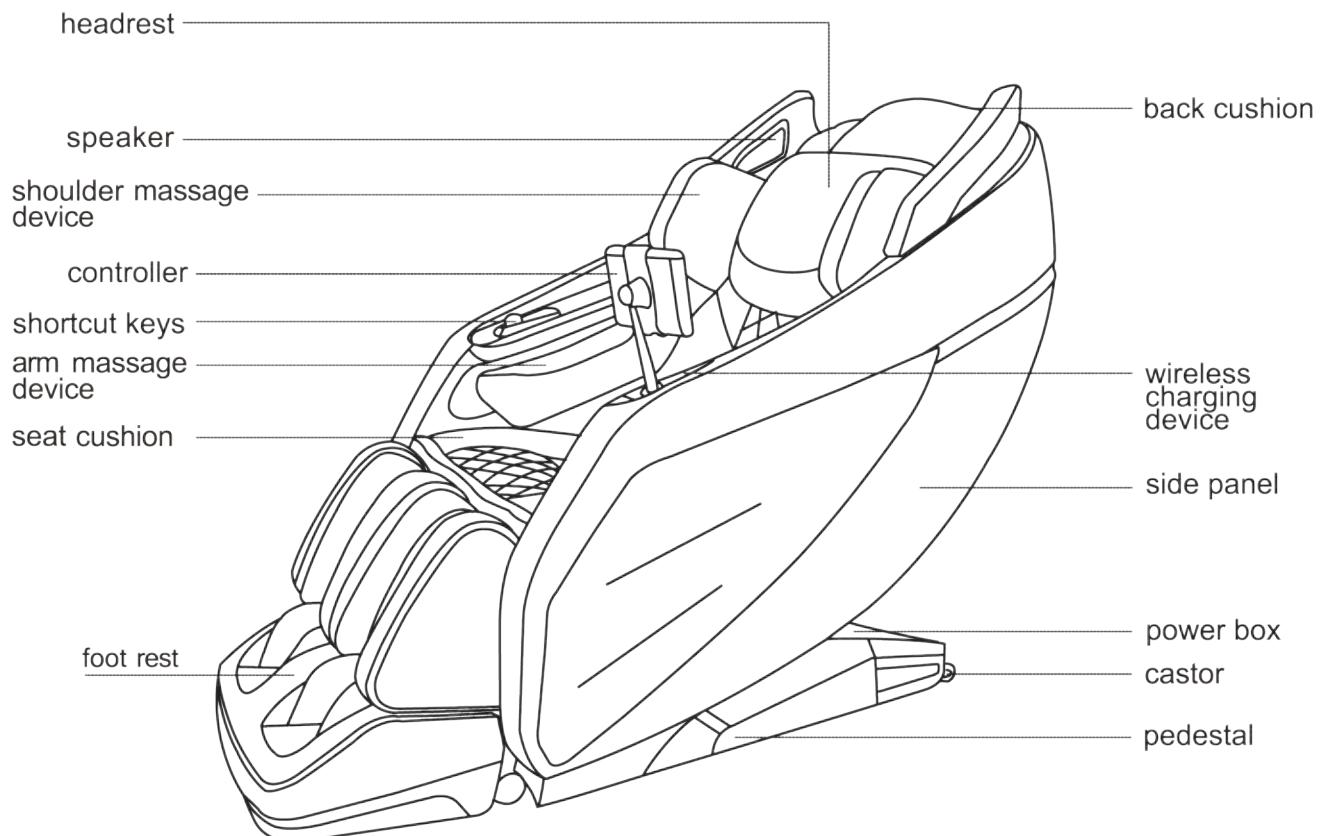
- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

NOTE 2: Any changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Refer to CC Part 18.213(c) and KDB 784748 Annex A section A.5:

- (a) Max. Electric Field Emissions: 6.14V/m
Max. Magnetic Field Emissions: 0.54 A/m
- (b) The maintenance of the system is the responsibility of the manufacturer.
- (c) To prevent interference, the product should be as far away from devices operating at the same frequency as possible.

FUNCTION OF EACH PART



Warning: Only rated supply 120V must be used to avoid damage to electronic elements.

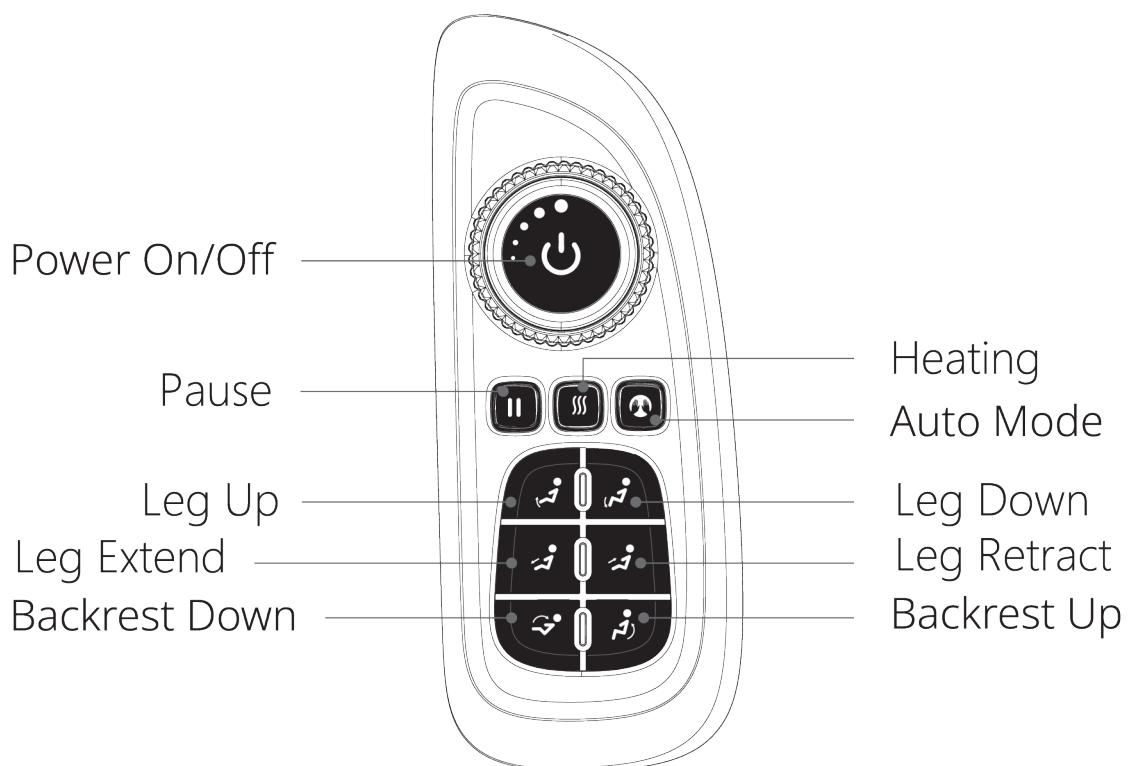
OPERATIONS

CONTROLLER MENU

Always press the Power Button to start massage.



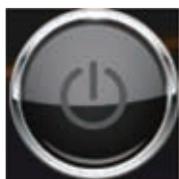
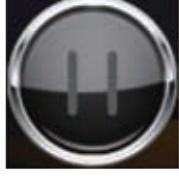
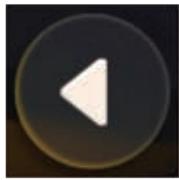
QUICK ACCESS CONTROL PANEL



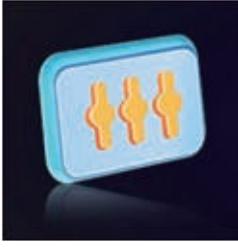
On/Off	Long press to turn the chair ON/OFF.
Pause	Long press to pause massage.
Heating	Long press to turn on the heating function.
Auto Mode	Short press to choose an Auto Mode. Every time this button is pressed, the auto massage mode will change in the sequence of Work Relief > Recovery > Deep Tissue > Night Care > Thai Stretch > Gentle > Athletic > Full Back > Circulation > Vitality > Neck & Shoulders > Waist & Hip
Backrest Up/ BackrestDown	Long press to adjust backrest.
Leg extend/ Leg retract	Long press to leg extend or retract.
Leg Up/ Leg Down	Long press to adjust legrest.

TABLET CONTROLLER INSTRUCTIONS

1. BASIC KEYS

KEYS	DESCRIPTION	GUIDELINES
	Off	Touch to stop massage and shut off the screen.
	Home	Touch to return to the home page.
	PAUSE	Touch once to pause current massage and the timer.
	RESTART	Touch to resume massage session, as well as timer.
	BACK	Touch to return to the previous page.
	ADJUSTMENTS	Touch to ADJUSTMENTS page.
	Air Massager	Touch to AIR MASSAGER page.

2. HOME PAGE MENU

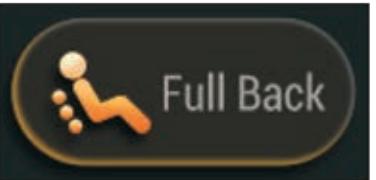
KEYS	DESCRIPTION	GUIDELINES
 AUTO	Auto	<p>Touch to choose from 12 AUTO massage programs.</p> <p>Work Relief , Recovery , Deep Tissue , Night Care , Thai Stretch , Gentle , Athletic , Full Back , Circulation , Vitality , Neck & Shoulders , Waist & Hip</p>
 MANUAL	Manual	<p>Touch to go to the MANUAL massage programs page, choose your preferred massage functions with desired intensities, fully customize your massage.</p>
 ADJUSTMENTS	Adjustments	<p>Touch to go to the ADJUSTMENTS page for zero gravity and massage intensities.</p>
 SETTINGS	Settings	<p>Touch to go to the system SETTINGS page for timer and other preferences.</p>

3. AUTO MASSAGE PROGRAMS



KEYS		GUIDELINES
	Work Relief	This is an exclusive massage program for those who have long sitting hours and business trips, as it can ease muscles and recuperate body.
	Recovery	A full body deep massage includes Shiatsu, Knocking and Tapping to ease and relax muscles after surgeries.
	Deep Tissue	A full body deep massage includes Shiatsu, Knocking and Tapping to ease and relax muscles after intense exercise, promoting the body circulation system to improve the muscle recovery after exercise.

	Night Care	<p>A program designed to provide a mild massage that is great for deep relaxation. This program helps improve the quality of sleep. The mild massage kneads body parts with strong intensity to light intensity to help soothe tired body gradually.</p> <p>Best used before going to sleep, the massage chair has the option of taking you into a Zer Gravity positionto to help you ease into complete relaxation.</p>
	Thai Stretch	<p>This massage will stretch your calf muscles and hamstrings with a combination of gentle tapping, kneading, rolling, and air compression.</p> <p>This massage is inspired by Thai Massage, which can effectively stretch legs and is the best choice to relax tired muscles and increase physical vigor.</p>
	Gentle	<p>A full body gentle massage includes Kneading and Heating, designed for daily use.</p>

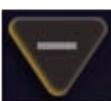
	Athletic	This auto program focuses on the tension points based on data collected from athletes, promoting the body circulation system to improve the muscle recovery after exercise.
	Full Back	A full body massage, mainly focusing on your full back with a deep air compression massage to relieve any back muscle tension.
	Circulation	A full body gentle massage includes Kneading and Heating, specially designed for women.
	Vitality	A full body strong massage includes Kneading, specially designed for men.
	Neck & Shoulders	A full body massage, that focuses on your neck and shoulders with a slow massage while incorporating soothing heat to relieve neck and back pain. The curved design of the backrest makes the massage rollers perfectly fit the shoulders and neck, effectively relieving fatigue and stiffness.

	Waist and Hips	A full body massage that mainly focuses on lower back, waist, and hips with a deep air compression massage to relieve full body pain.
---	----------------	---

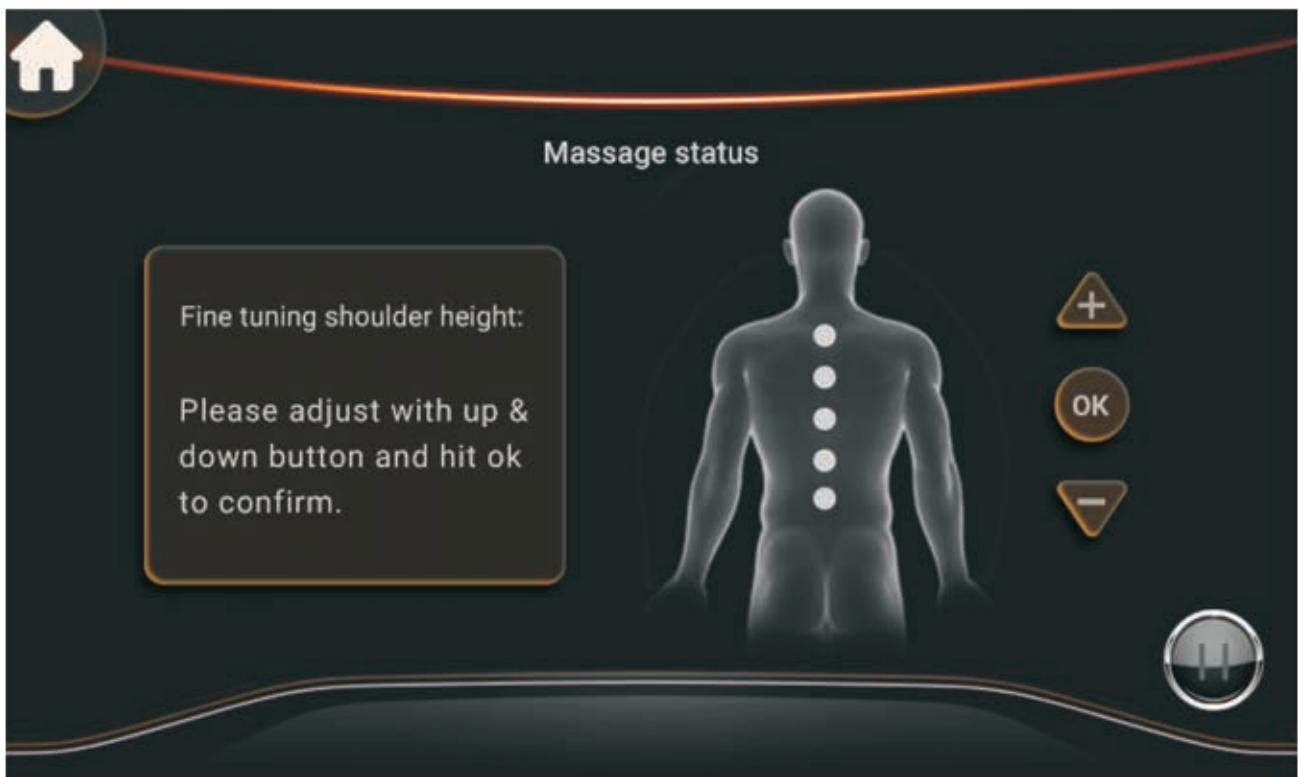
Backrest Scanning Instructions

A body scan is performed every time the power is on and the user selects an Auto Mode, it is designed to ensure an accurate and personalized massage for every individual.

After the power is on, sit on the chair slowly with your back against the chair back, and your head against the pillow so that the massage sensors can detect your Shoulder Height and the whole back.

The Tablet Control gives beeping sounds when completing a body scan. When micro-adjustment is needed, just touch Roller up/down  , then  to start massage. Make sure the body scan is finished completely for the best massage.

Once body scan is completed, it will skip to the massage status page without further actions.

**Note:**

Heating, foot rollers can be turned on/off manually during any auto mode, see below.

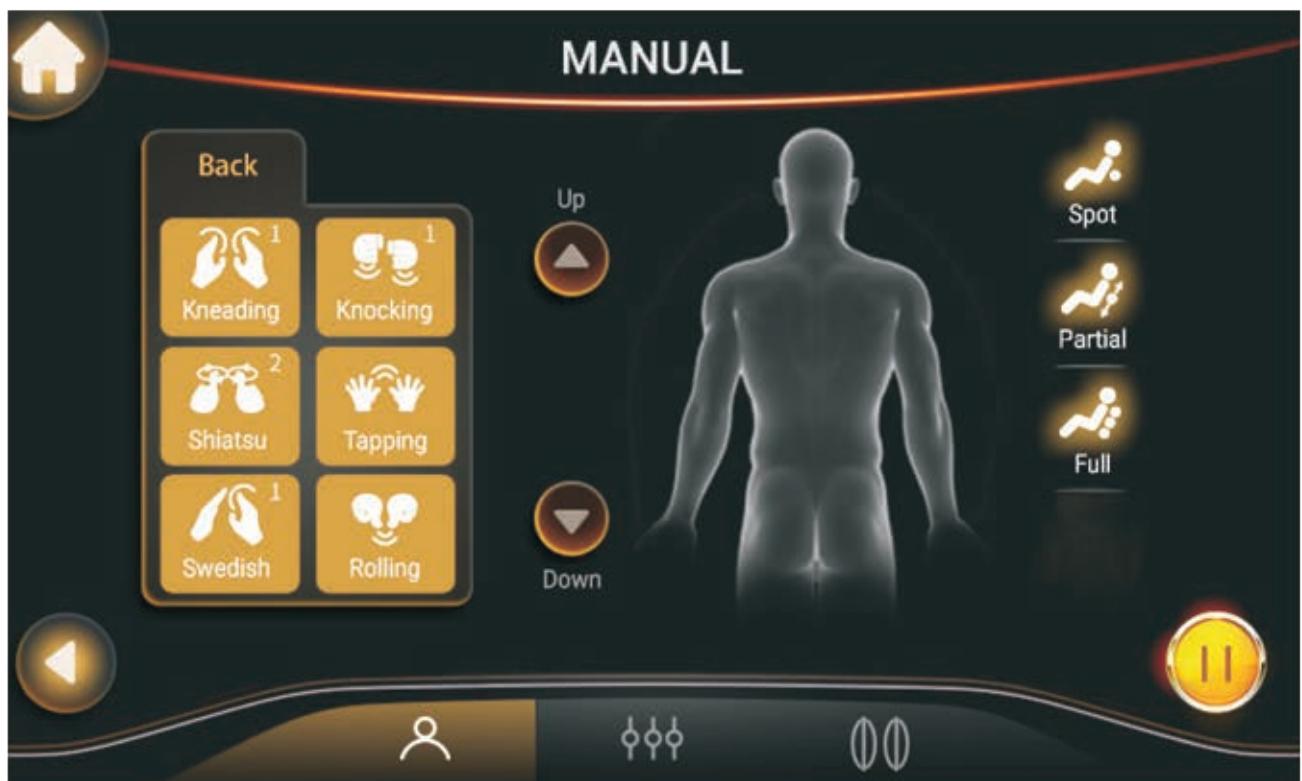


4. MANUAL MASSAGE PROGRAMS

6*2 massage functions on demand.

Kneading, Shiatsu, Swedish, Knocking, Tapping, Rolling.

- Touch to choose a manual massage function for a full back massage at medium strength.
- Manual setting modes can not be combined with other modes.
- The default setting is spot-specific massage with medium width.

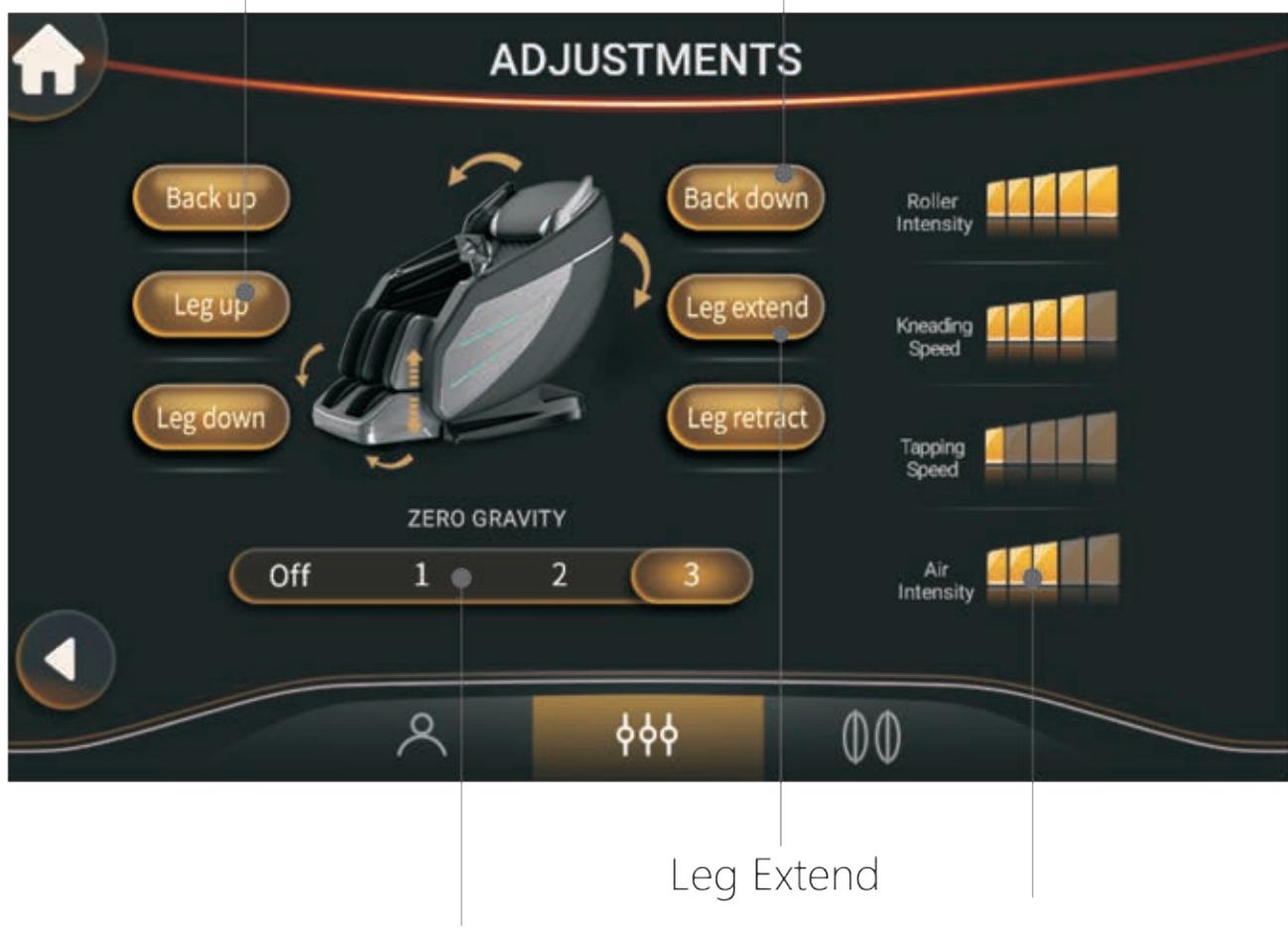


KEYS	DESCRIPTION	GUIDELINES
 	Roller Spot/ Partial Massage	<p>Touch Spot to activate localized Spot Massage while the original massage mode continues to run in the current position.</p> <p>Touch Partial to activate Partial Massage, rollers take the latest massage spot position as the center, and will travel approximately 3 inches up or down. Only effective under manual massage mode.</p>
	Full Back Roller	Full Back Massage.
 	Roller Up or Down	Up or down adjustment.

5. ADJUSTMENTS

Legrest Angle Adjustment

Backrest Angle Adjustment



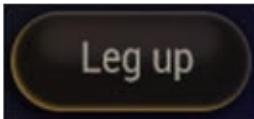
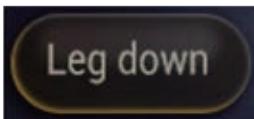
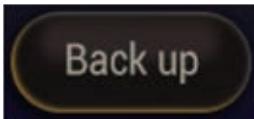
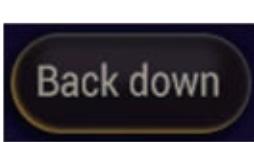
Leg Extend

Zero Gravity Angle Adjustment

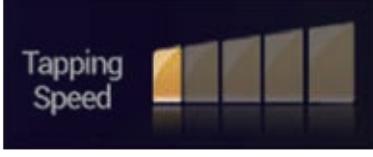
Intensity Adjustment

Zero-gravity is based on technology originally pioneered by NASA. Using the back, seat, and the leg rest, the chair puts the heart and knee at the same level, lowering blood pressure. Not only that, but all the gravitational forces can be dispersed to the massage chair, reducing heart and spine pressure.

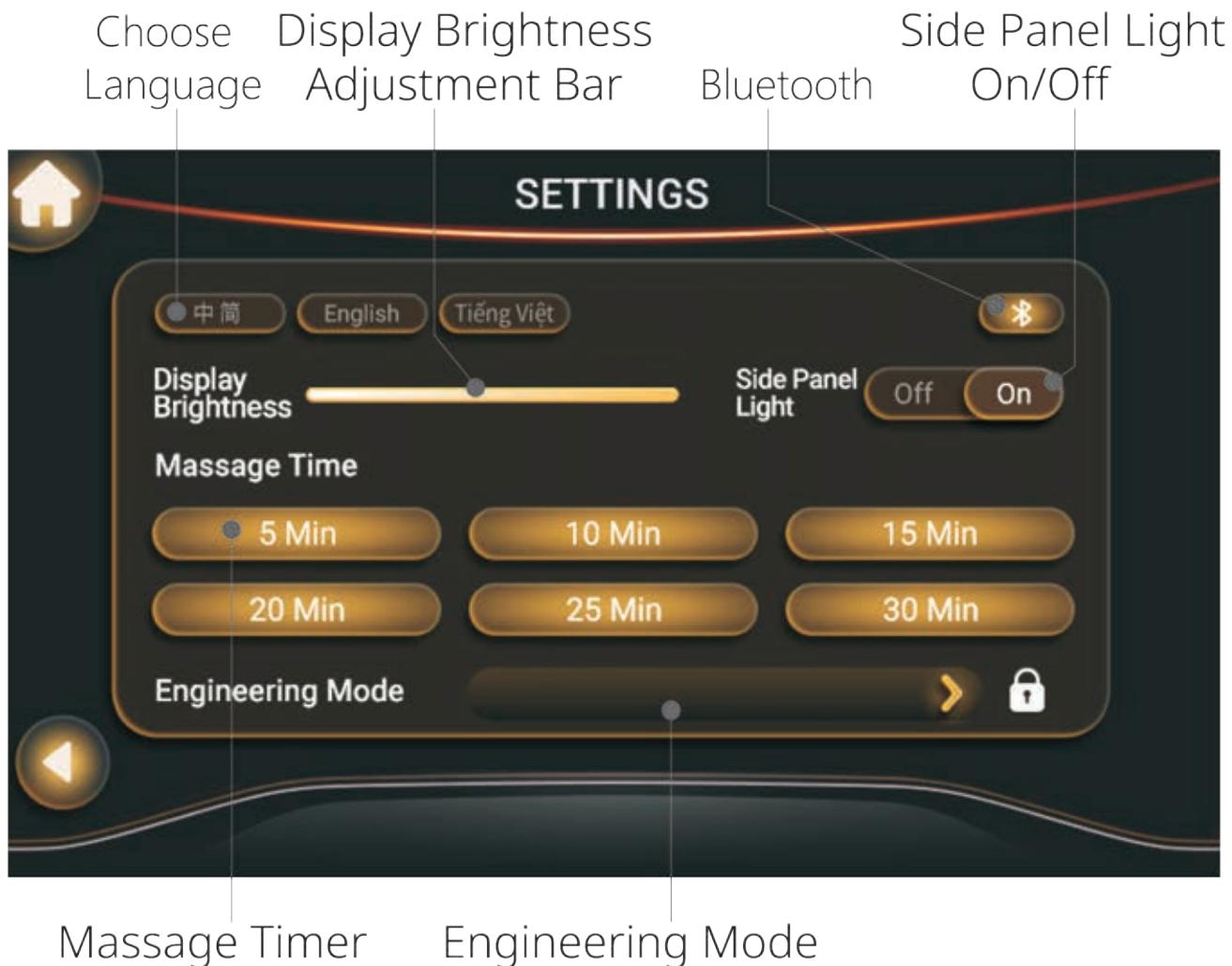
Zero Gravity Angle Adjustment

KEYS	DESCRIPTION	GUIDELINES
 	Angle Adjustment Key for Legrest	<p>Leg Up Long-press this button, to raise the legrest. Let go of the button to stop raising the legrest.</p> <p>Leg Down Long-press this button, to lower the legrest. Let go of the button to stop lowering the legrest.</p>
 	Angle Adjustment Button for Backrest	<p>Back Up Long-press this button, to raise the backrest. Let go of the button to stop raising the backrest.</p> <p>Back down Long-press this button, to lower the backrest. Let go of the button to stop lowering the backrest.</p>
 	Leg extend or retract	<p>Extend Legrest Long-press to extend legrest.</p> <p>Retract Legrest Long-press to retract legrest.</p>
	Zero Gravity Function Control	Touch this button to activate the zero gravity position, and choose from Zero Gravity mode 1, Zero Gravity mode 2 and Zero Gravity mode 3.

Other Adjustments

KEYS	DESCRIPTION	GUIDELINES
	Roller Intensity	SLOW / MEDIUM / FAST Effective only under manual massage mode.   
	Kneading Speed	5 levels kneading speed adjustments.
	Tapping Speed	5 levels tapping speed adjustments.
	Air Intensity	5 levels tapping air intensity adjustments.

6. SETTINGS



BLUETOOTH INSTRUCTIONS

Play Mobile Music/Bluetooth speaker:

Turn on Bluetooth of mobile, then select Bluetooth device of massage chair (OSAXXXX) Once connected, any audio may be played through the speakers using the mobile device.(Note: no password needed)

SAFETY INSTRUCTIONS

1. Should any abnormality occur, shut off and detach power plug from socket.
2. When power cord damaged, it must be changed, please purchase only from accredited reseller.



Warning: Do not dismantle, modify or change components to avoid electric shock or fire.

DAILY CARE AND MAINTENANCE

Storage guide

- © Do not dismantle back of the chair. Electrical shock may happen when in touch with the inner elements of chair.
- © An unauthorized person must not dismantle the chair to avoid electric shock or fire.
- © Make sure to shut off and detach the power plug when cleaning to avoid electric shock or hazard.
- © Keep the massage chair clean and free of dust. Do not place the massage chair in a damp environment.
- © If not in use for a long period, please cover it with a cloth to free dust.
- © Do not place the massage chair under sunlight.

Maintenance

- © Use neutral detergents only
- © Use neutral detergents to clean the back cushion, head cushion, and leg cushion. Let it dry.
- © Use detergent and cloth to clean the plastic, pipe and leather material of massage chair.
- © Make sure not to use organic solvents as gasoline, amyl acetate, benzene for cleaning to avoid fading of the color, aging, or breakage.
- © Do not spray pesticide to avoid electric shock or fire.
- © Please clean with dry cloth to avoid mechanical failure.

Controller maintenance

- © Do not place the massage chair near heat or fire.
- © Keep away from dampness to avoid electric shock.
- © Please clean dust with a soft cloth or a small banister brush.
- © Make sure not to use violent detergent such as benzene, gasoline, amyl acetate.

TROUBLE SHOOTING

Checking is necessary for frequent usage.

Make sure the machine works properly before usage, after not in use for a long time.

Abnormality situation	Solution guide
Failure to start and indicator light is off	<ul style="list-style-type: none">© Check that the plug is well connected© Check if the switch on the chair is turned to "I" position
Sound in operation	<p>The sound is normal in following situation</p> <ul style="list-style-type: none">© Sound for kneading or knocking© Sound changes after speed adjustment© Sound from the box in massage© Vibration sound for backrest in massage© Sound for air input and output© Sound for airbag inflating and venting© Sound for leg adjustment© Sound for motor when the massage chair is loaded

SPECIFICATION

Product name	Massage Chair
Model	OS-3D Champ 2
Power supply	120V~ 60Hz
Rated power	150W
Rated working time	15 min